L I F E I N B A L A N C E

SOEHNLE

MASSAGE CUSHION

CHILL & GO 3D

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OPERATING INSTRUCTIONS

MODE D' EMPLOI

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BRUKSANVISNING

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Händler

IMPORTANT NOTES!

PLEASE STORE THESE INSTRUCTIONS FOR FURTHER USE!

1. Proper use	13
2. Safety information	14
3. Care instructions	16
4. Scope of delivery	16
5. Product description	17
6. Start-up Operation	18
7. Utilization	19
8. Problems – Causes - Troubleshooting	20
9. Storage, transport, and periodic inspection	20
10. Technical Specifications	21
11. Declaration of conformity	21
12. Disposal	21
13. Hotline	21
14. Guarantee	22

1. Proper use

We thank you for your purchase of this device.

With this portable massage seat pad you will be able to have a massage at anytime and anywhere. The massage seat pad can be used on an office chair or any other seat or chair with a backrest of sufficient height and a sufficient seat cushion depth.

Massage is one of the oldest healing methods known to mankind. Massages usually serve to relax muscles and increase your feeling of well-being. Massages with this device are relaxing or vitalizing and have proven effects on muscle tensions and fatigues.

The device is not a replacement for medical or physiotherapeutic treatment. Please consult your physician if you are unsure whether this massager is suitable for you.

The device is intended exclusively for domestic use, not for commercial purposes, the use in hospitals or doctor's offices.

In order to be able to utilize all the benefits of your massager, please thoroughly read these instructions and observe the notes on care and safety. You will then be able to enjoy this massager for a long time.

Please store these instructions for further use or in case that you hand this device over to another person.

Further information can be found at our website at www.soehnle.com.

2. Safety information

Danger

Improper use may result in hazards that may lead to death or severe injury.

Warning

Improper use may, in certain circumstances, cause injuries and pose risks to your health.

Attention

Improper use may result in injury or property damage.

Notes

Important information

Danger

- This device uses mains voltage. Therefore, the same safety regulations apply as for any other electric device:
 - Operate the device only with the stated supply voltage
 - Inspect all components of the device for faultless condition prior to each start-up of the device. In the event of a defect or malfunction, switch the device off immediately and disconnect the power plug from the power network. Repairs must be performed exclusively by authorized customer service technicians. Never open the device. The individual components of this device may not be disassembled, repaired or modified.
- Use the device only in dry rooms and with a dry body – never in the bathtub, shower, above filled sinks, in pools, saunas, or outdoors. 
- Never operate the device in the close vicinity of heat sources.
- Do not pull the power cord across sharp edges.
- Do not pull on the cord or the device in order to disconnect the power plug or to transport the device.
- Disconnect the power plug if the device is not in use.
- If you dispose of the device at any one point in time, please cut off the cable to make it useless.
- This device may not be used by persons (to include children) with physical, sensory, or mental disabilities or by persons with insufficient experience or knowledge unless such persons

use the device under the supervision of a responsible person or after receiving necessary instructions.

- Please ensure that children do not use the device as a toy.
- Disconnect the device from the power network prior to cleaning it. When cleaning, please ensure that no liquids enter the switches and buttons and that the device is used only in an absolutely dry condition.
- Do not disassemble the device. It does not contain any parts that can be repaired by the user.
- The cover of the device, i. e. the zipper, cannot be opened; this means that a repair of interior components must be performed by specialized technicians.
- Do not insert any needles or pointy objects into the massager
- Never use this device with other electronic medical devices, electronic life-support systems, e. g. heart/lung machines, or mobile electronic medical devices like electrocardiographs. This device may cause malfunctions in the above-mentioned electronic medical devices.

Warnings

- Consult a physician prior to use if you
 - are pregnant
 - have a pacemaker or similar
 - suffer from an illness (e. g. flu, fever, varicose veins, thrombosis, jaundice, diabetes, infections, sciatica, herniated discs, or similar)
 - are recovering from surgery.
- Never use the device with open wounds, burns, hematoma, swellings, or similar conditions.
- Never treat joints, the spine, head, and other sensitive body parts.
- Never use the device under the influence of alcohol.
- Stop the massage if it feels painful or uncomfortable to you.
- If you feel pain in muscles or joints after using the device for a long period of time, immediately stop using the device and consult a physician.
- For the same reason, the device should never be used underneath pillows, blankets, or similar. Never use the device in a bed to prevent falling asleep while the device is in use.

- To prevent any hazards if the power cord is damaged, it must be replaced by the manufacturer, its customer service or a similarly qualified person.
- The device is equipped with heated areas. Persons that are insensitive to heat must be careful when using the device.
- The device is equipped with heated areas. Persons that are insensitive to heat must limit the use to 15 minutes or may use the device only under supervision.
- Never use the device without a cover.

Attention

- Ensure that all components can move freely during operation.
- Use the device only in a vertical position, never while laying down, e. g. on a bed or couch, as this will cause excessive pressure on the body.
- Never use while driving a vehicle.

Notes

- If the device is used just prior to going to bed, please remember that a massage may be revitalizing and cause sleep disorders.
- Automatic off after 15 min. This corresponds to the optimum duration of a massage session that has a revitalizing effect on the body.

3. Care instructions

- Only use a cloth moistened with mild detergent for cleaning.

Danger

- Disconnect the power plug from the power outlet and let the device cool off prior to cleaning the device.
- Only use the device after it has dried completely.
- The cover of the device, i.e. the zipper, cannot be opened; this means that cleaning of interior components must be performed by specialized technicians.

**Attention**

- Never use abrasive and corrosive cleaning agents and never submerge the device in water.

**4. Scope of delivery**

Massage seat pad with seat pad, power supply unit, and operating instructions

5. Product description

Operating element


Functional keys

1. On/Off switch
2. Heat function key

Functional keys for

3D Shiatsu deep massage (left)

3. Full back
4. Upper back
5. Lower back
6. Selective massage toward the top
7. Selective massage toward the bottom

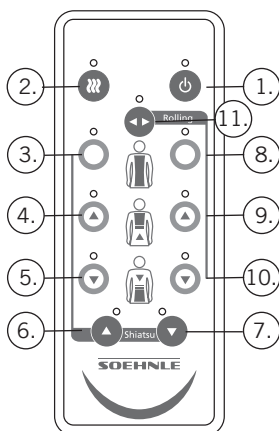
Functional keys for

Rolling massage (right)

8. Full back
9. Upper back
10. Lower back
11. Width regulation

Massage seat pad

12. Massage seat pad
13. Removable backrest pad with hook-and-loop fastener
14. Pocket for operating element
15. Attachment strap with hook-and-loop fastener
16. Operating element
17. Power cord with power supply unit
18. 4 rotating massage heads
19. Removable seat pad with hook-and-loop fastener
20. Carrying handle

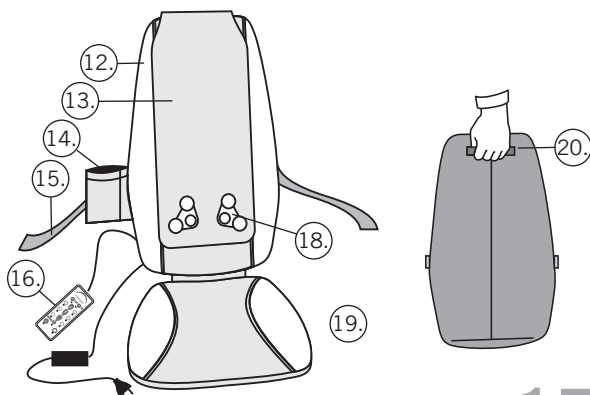6. Start-up Operation

Please read through these operating instructions in proper order.

Attaching the massage seat pad

The massage seat pad can be used on an office chair or any other seat or chair with a vertical backrest of sufficient height and a sufficient seat cushion depth.


The convenient hook-and-loop fasteners will ensure that the massage seat pad does not slide during use.

Connecting the power supply

Connect the power plug of the power supply unit to the power outlet and the low-voltage plug to the low-voltage socket on the massage seat pad.

Note

For safety reasons the massage seat pad may be operated only with the enclosed power supply unit, not any other device.

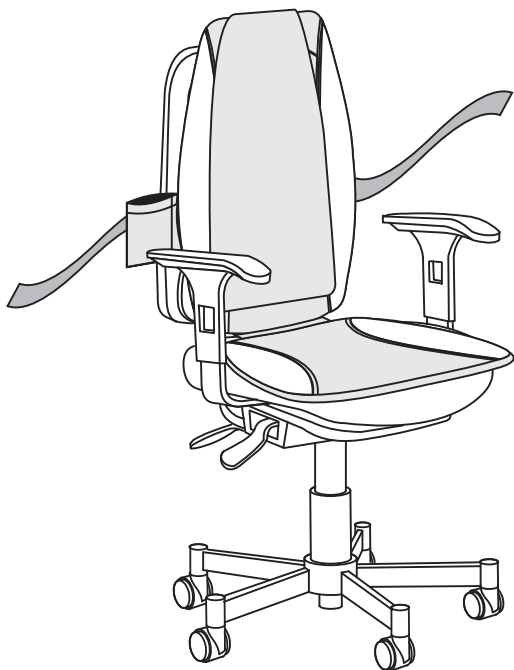7. Utilization

Test the massage seat pad prior to use by carefully sitting down on it. When in resting position, the massage heads will be in your lumbar area. Carefully lean onto the back rest to find out how much pressure you would like to use.

It is important that you sit in the middle so that the massage heads do not massage the spine. Due to the removable seat pad, the massage seat pad can easily be positioned on all commonly available chairs.

On/off switch

Use the On/Off switch to switch the device on.

Selecting the massage type and massage area

The device is equipped with 2 different massage types, each with a variety of options.

1. 3D deep tissue shiatsu massage

"Shiatsu" is a Japanese term meaning "finger pressure": shi = finger, and atsu = pressure. It is a vitalizing and strengthening, circular finger pressure massage.

2. Rolling massage

A rolling massage is a calming and relaxing massage along the spine.

Use the operating keys of the desired massage type to select the area to be massaged.

You can then choose between a massage of the full back, the upper or lower back.

The massage type "Shiatsu" can also provide a selective massage. The precise point to be massaged can be controlled by pushing and holding the ▲/▼ keys.

During the "rolling" massage it is also possible to regulate the individual width of the massage heads by pushing and holding the ◀ ▶ key.

After selecting the massage area, you will be able to add pleasant warmth by using the heat function key.

Note: The heat function can be used only in combination with a massage type.

The massage strength can be individually adjusted with the removable backrest.

The recommended duration of a back massage is 15 minutes. Therefore, the device will automatically switch off after that time.

Switching off

The device will automatically switch off after pushing the On/Off key or after 15 minutes.

The massage heads will return to the starting position.

Disconnect the power plug from the power network.

8. Problems – Causes - Troubleshooting

Problem	Cause	Troubleshooting
Massage heads rotate only slowly	Excess load on massage heads	Use less force when pressing the body into the backrest
Device cannot be switched on	The device is not connected to the power network.	Connect power plug and switch device on
Device cannot be switched on, even though indicator light on power supply unit lights up	The device might not be switched on	Use the On/Off switch to switch the device on and select a massage area.
	The overload protection may have been triggered.	Disconnect power plug and reconnect after a few minutes

9. Storage, transport, and periodic inspection

- Store the massager in a dry location if it is not in use.
- Let the device cool off prior to storing it.
- Regularly inspect the device for possible signs of wear or damage.
- Please return the device to your point of sale in case of any wear or malfunction. Your retailer will assist you.
- Do not place any heavy objects on top of the device while storing it in order to prevent possible damage to components. Do not stand on the device.
- Do not drop the device. This may result in damages.
- Use care when transporting the massager.

10. Technical Specifications

Dimensions	~ 26.2 (L) / 17.7 (W) / 2 (H) inch
Weight	~ 9.24 lb
Ratings	for adapter: input: 100 - 240 V, 50/60 Hz, 1.5A output: 12V === 3.0A for massage: input: 12V === 3.0A Total input power (adapter + massage): 50 W

The device is equipped with double protective insulation and corresponds to protection class II.

11. Declaration of conformity

SOEHNLE herewith declares that this device is in compliance with the basic requirements and all other applicable stipulations of Directives 2004/108/EC and 2006/95/EC.

The EU Declaration of Conformity can be found at www.soehnle.com.

12. Disposal of used electronic devices

After its service life, please dispose of this device according to the Directive on Waste Electrical and Electronic Equipment 2002/96/EC (WEEE).

The symbol on the product or its packaging informs you that this product cannot be treated as normal household waste, but must be returned at a collection point for recycling electronic and electrical devices.

13. Customer service

Our customer service will be happy to assist you if you have any questions or comments:

Germany Ph.: (08 00) 5 34 34 34 (free of charge)
International Ph.: +49 26 04 97 70
Monday through Friday 08:30 a.m. to 12:00 noon CET

If you'd like to contact us in writing, please use the contact form at the SOEHNLE homepage www.soehnle.com or use postal mail.

14. Guarantee

SOEHNLE guarantees rectification of faults based on material or production errors, free of charge, either by repair or replacement, for 3 years as of the date of purchase. Please retain the receipt of purchase and guarantee slip. In case of a guarantee claim, please return the device, guarantee slip and receipt of purchase to your seller.

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Leifheit AG